

1.16


A Typical Workflow of Routine Analysis

Std. Familiarization Time: 0:10:00

Approximate Video Time: 0:05:00

[Play the Video](#)

MassHunter Basics for Agilent ICP-MS
Rev.3.0.2 – February 2021
71

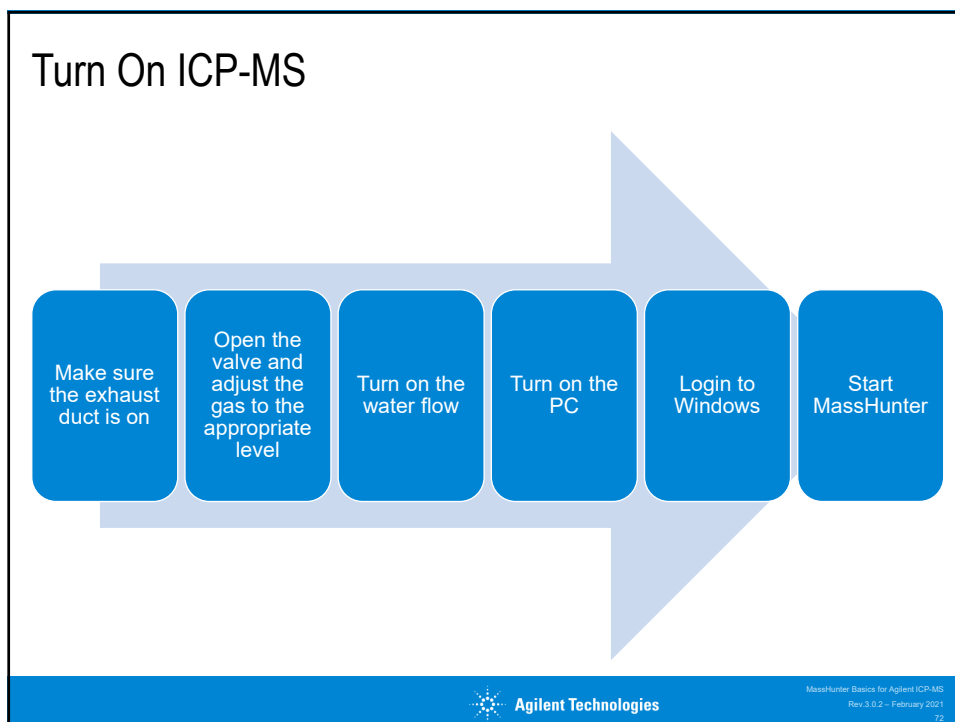

Plasma On

Purge the
cell gas
line

Preparing
required
solution

Install the
peripump
tubes and
adjust the
tension

Set up
Startup

Move the
probe into
the DI
water

Turn on
the
plasma

Running a Batch

Create a
batch from
a batch
template

Perform
an auto
tune

Set up the
sample list

Check the
on-line
ISTD
signal

Add the
batch to
run

Check the
analysis
results

